




Aripiprazole Oral Solution – BP 2026

These chromatograms are provided for information only as an aid to analysts and are intended as guidance for the interpretation and application of BP monographs.


Typical spectrum for the Identification test for Aripiprazole Oral Solution as published in BP 2026.

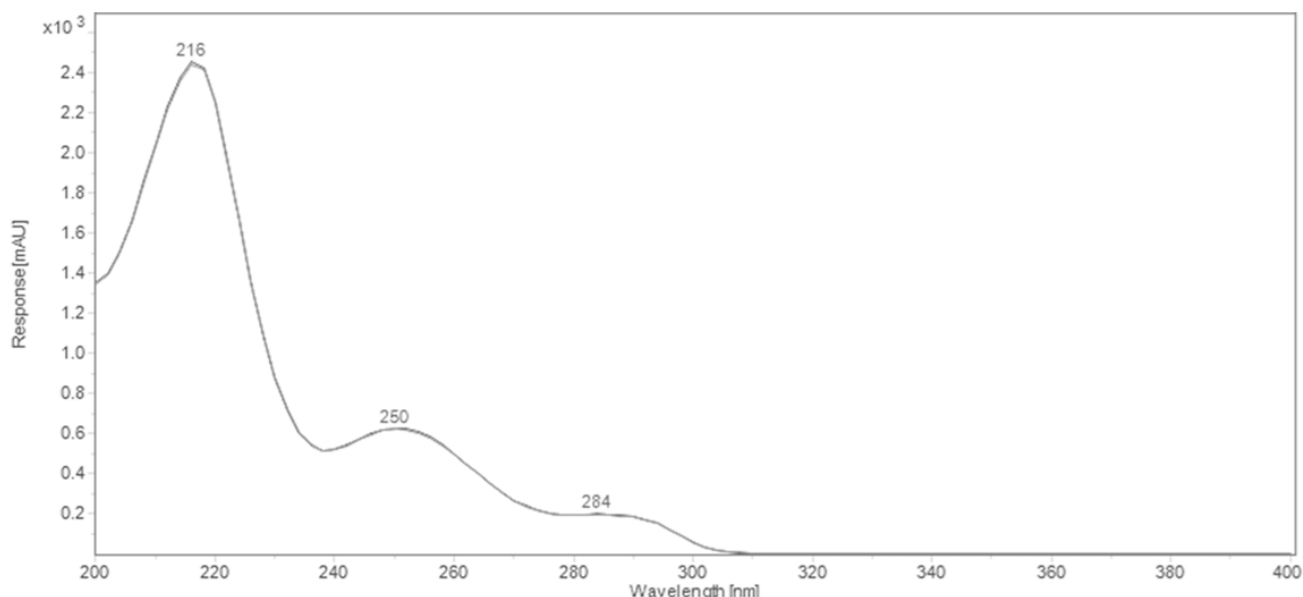

Column	Poroshell EC-CN (100 mm x 3.0 mm, 2.7 µm)
Method Ref.	Identification for the Aripiprazole Oral Solution monograph from BP 2026
Buffer	0.264% w/v diammonium hydrogen orthophosphate in water, adjusted to pH 3.0 with orthophosphoric acid
Mobile Phase	Acetonitrile: buffer (30:70, v/v)
Diluent	Dichloromethane: acetonitrile: 0.2% v/v orthophosphoric acid (0.05:0.95:9, v/v/v)
Flow rate	0.6 mL/min
Column Temp	25°C

Injection Volume	2 µL
Detection	230 nm

Typical chromatogram for solution (3) from the Related Substances test for Aripiprazole Oral Solution as published in BP 2026.

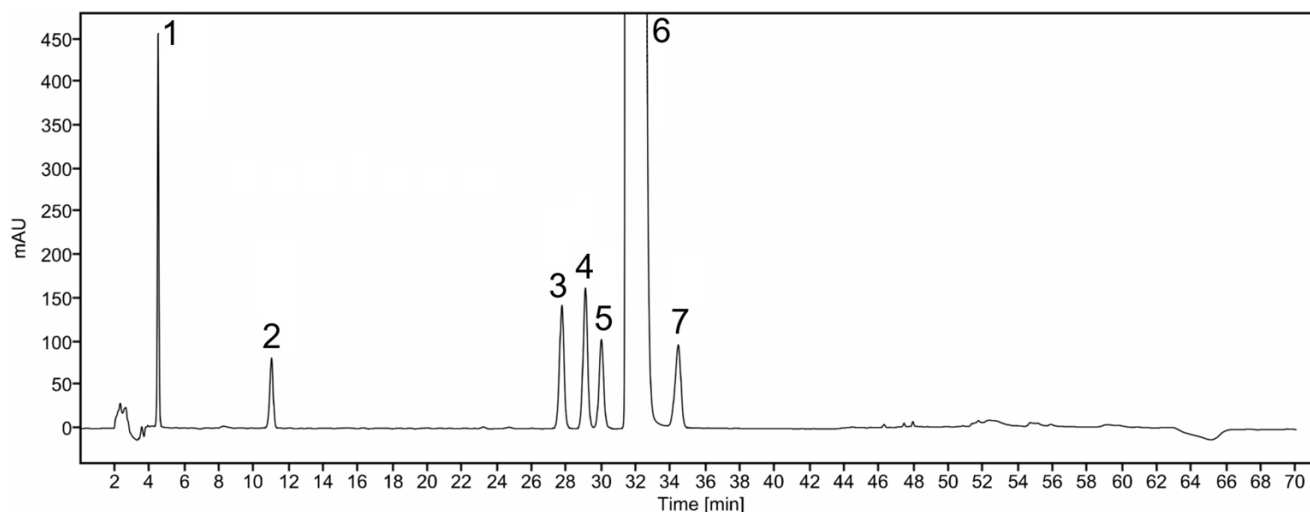


Peak ID: 1: Impurity A. 2: Impurity B. 3: Impurity E. 4: Impurity 2. 5: Impurity 1. 6: Aripiprazole. 7: Impurity F.

Column	Luna C18 (2) (250 mm x 4.6 mm, 5 µm)
Method Ref.	Related Substances for the Aripiprazole Oral Solution monograph from BP 2026
Mobile Phase A	0.111% w/v potassium dihydrogen orthophosphate and 0.12% w/v sodium pentanesulfonate, pH adjusted to 3.0 with 10% orthophosphoric acid
Mobile Phase B	Acetonitrile
Solution A	Equal volumes of acetonitrile and methanol
Solution B	0.111% w/v potassium dihydrogen orthophosphate, previously adjusted to pH 3.0 with 10% v/v orthophosphoric acid
Diluent	Solution A: 0.1M hydrochloric acid: methanol: solution B (40:1.5:6:2.5, v/v/v/v)
Flow rate	1.0 mL/min
Column Temp	25°C
Injection Volume	20 µL
Detection	215 nm
Gradient	

Time (minutes)	Mobile phase A (% v/v)	Mobile phase B (% v/v)	Flow rate (mL/min)	Comment
0 – 10	72 → 70	28 → 30	1.0	linear gradient
10 – 40	70 → 60	30 → 40	1.0	linear gradient
40 – 50	60 → 37	40 → 63	1.0	linear gradient
50 – 60	37	63	1.0	isocratic
60 – 61	37 → 72	63 → 28	1.0	linear gradient
61 - 70	72	28	1.0	re-equilibration

Typical chromatogram for solution (2) from the Assay test for Aripiprazole Oral Solution as published in BP 2026.

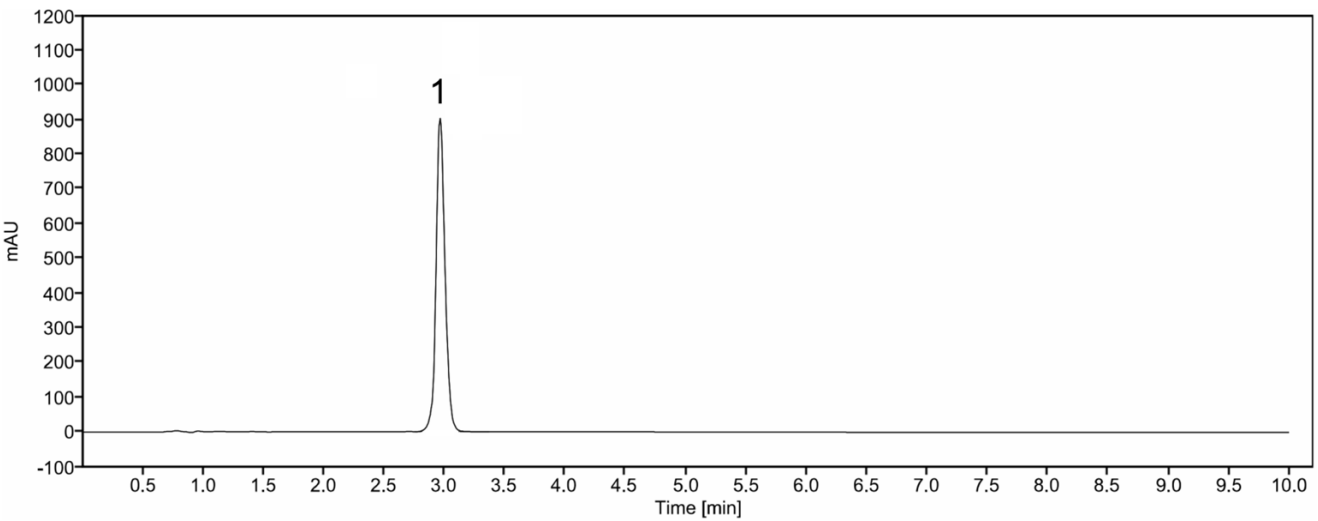

Peak ID: 1: Aripiprazole

Column	Poroshell EC-CN (100 mm x 3.0 mm, 2.7 µm)
Method Ref.	Assay for the Aripiprazole Oral Solution monograph from BP 2026
Buffer	0.264% w/v diammonium hydrogen orthophosphate in water, adjusted to pH 3.0 with orthophosphoric acid
Mobile Phase	Acetonitrile: buffer (30:70, v/v)
Diluent	Dichloromethane: acetonitrile: 0.2% v/v orthophosphoric acid (0.05:0.95:9, v/v/v)
Flow rate	0.6 mL/min
Column Temp	25°C
Injection Volume	2 µL
Detection	230 nm