




Nevirapine Prolonged-Release Tablets


Typical chromatogram for solution (2) in the Assay test for Nevirapine Prolonged-Release Tablets as published in BP 2023.

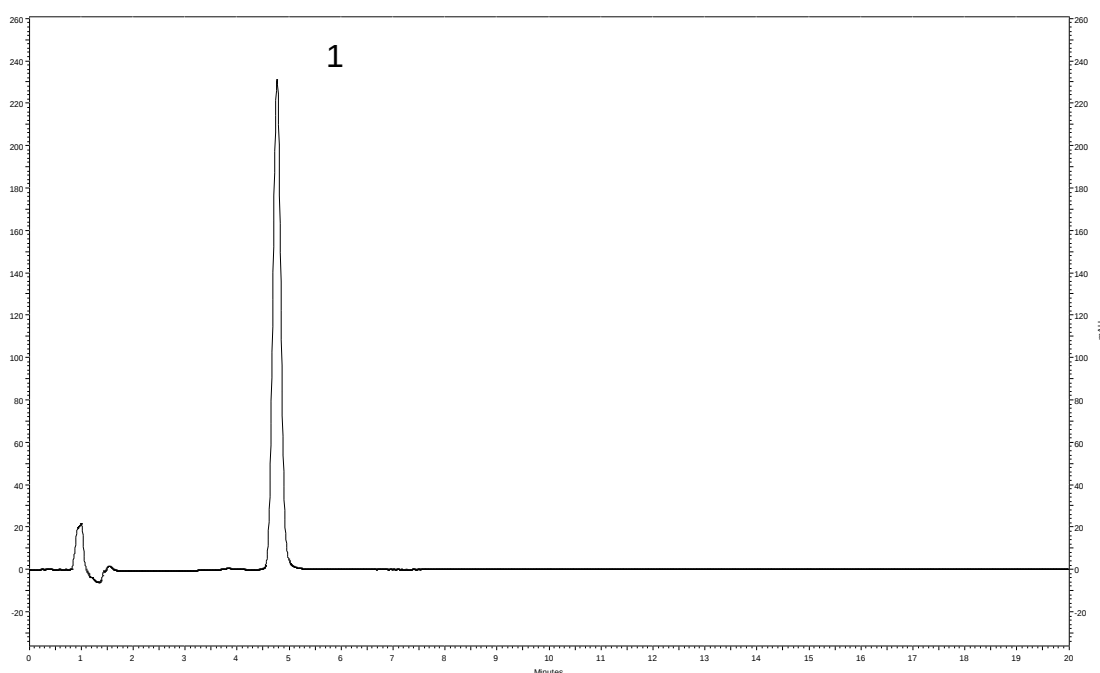

Peak ID: 1: Nevirapine

Column	Waters, Symmetry C18 (150 mm x 3.9 mm, 5 µm)
Method Ref.	Assay for the Nevirapine Prolonged-Release Tablets monograph from BP 2023
Mobile Phase	Acetonitrile : water (23:77, v/v)
Diluent	Methanol : water (80:20, v/v)
Flow rate	1.0 mL/min
Column Temp	25 °C
Injection Volume	20 µL
Detection	214 nm

Typical spectrum for solution (2) in the Identification B test for Nevirapine Prolonged-Release Tablets as published in BP 2023.

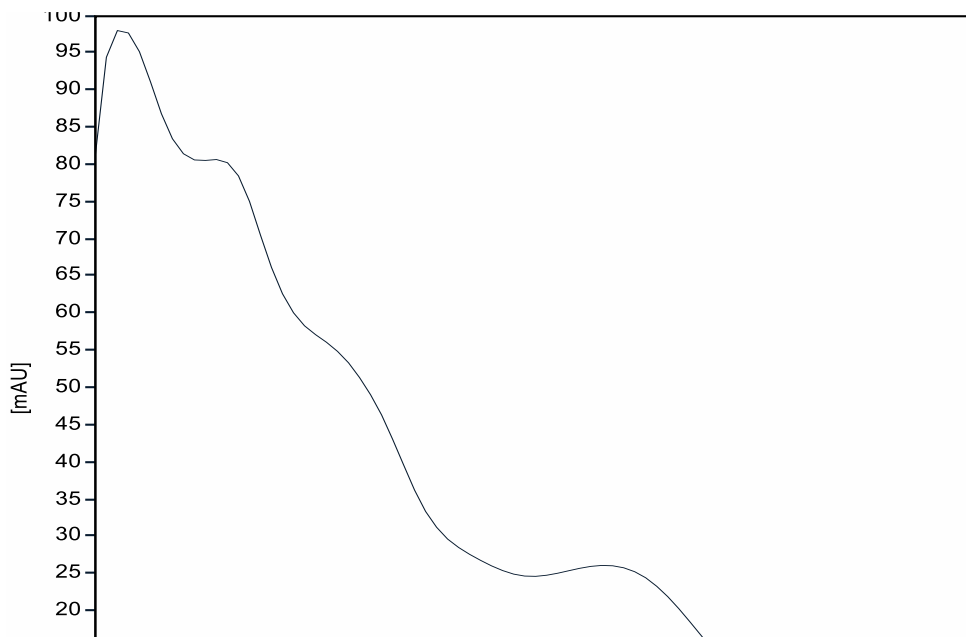

Column	Waters, Symmetry C18 (150 mm x 3.9 mm, 5 µm)
Method Ref.	Assay for the Nevirapine Prolonged-Release Tablets monograph from BP 2023
Mobile Phase	Acetonitrile : water (23:77, v/v)
Diluent	Methanol : water (80:20, v/v)
Flow rate	1.0 mL/min
Column Temp	25 °C
Injection Volume	20 µL