




Warfarin Tablets – BP 2027

These chromatograms are provided for information only as an aid to analysts and are intended as guidance for the interpretation and application of BP monographs.

Typical chromatogram for solution (3) from the Related Substances test for Warfarin Tablets as published in BP 2027.

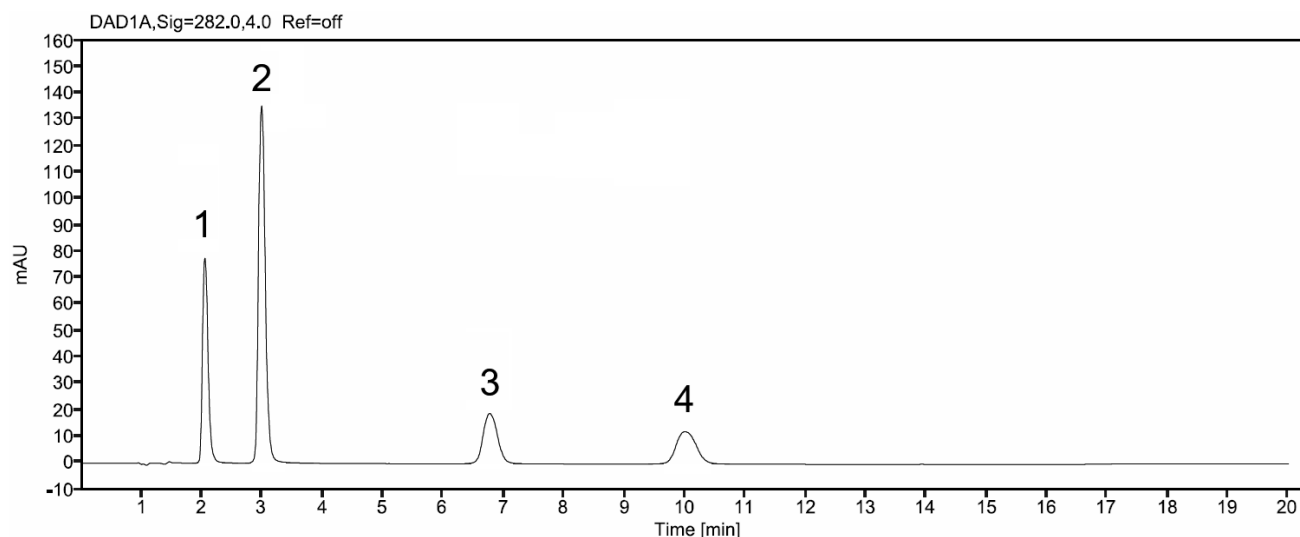

Peak ID: 1: Impurity B. 2: Impurity C. 3: Warfarin. 4: Impurity A.

Column	Zorbax Phenyl-Hexyl (150 mm x 4.6 mm, 3.5 µm)
Method Ref.	Related Substances for the Warfarin Tablets monograph from BP 2027
Mobile Phase	1 volume of glacial acetic acid, 1 volume of tetrahydrofuran, 42 volumes of water and 57 volumes of methanol
Diluent	50% v/v methanol
Flow rate	1.5 mL/min
Column Temp	35°C
Injection Volume	100 µL

Detection	282 nm
------------------	--------