




Spironolactone Oral Suspension – BP 2016


These chromatograms are provided for information only as an aid to analysts and intended as guidance for the interpretation and application of BP monographs.

Typical chromatogram for solution (4) in the Related Substances test for Spironolactone Oral Suspension as published in BP 2016.

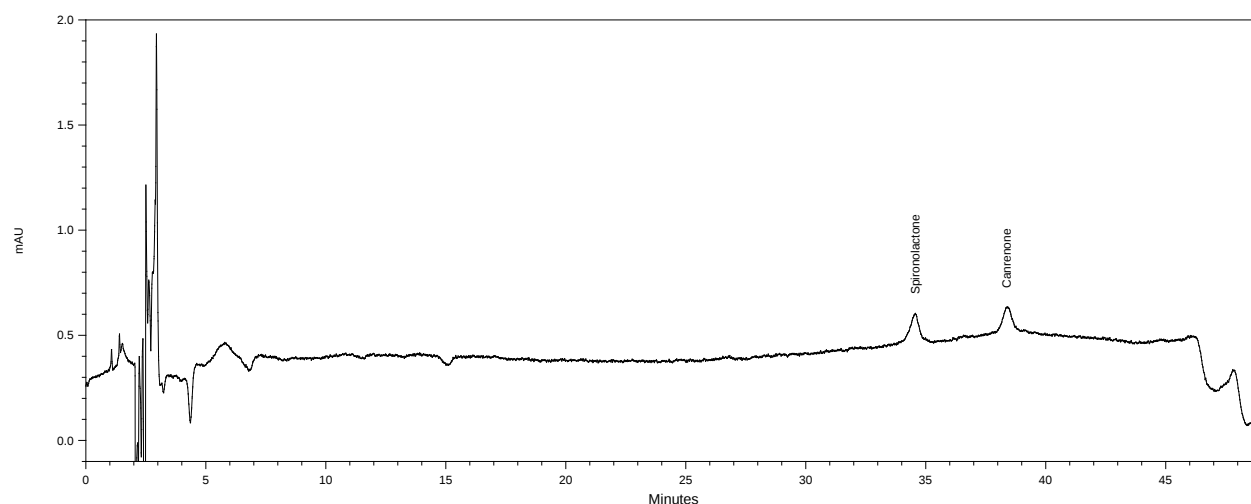

The chromatogram shows two peaks spironolactone (>1 mAU) at 34.6 minutes and canrenone (>1 mAU) at 38.4 minutes.

Typical chromatogram for solution (1) in the Related Substances test for Spironolactone Oral Suspension as published in BP 2016.

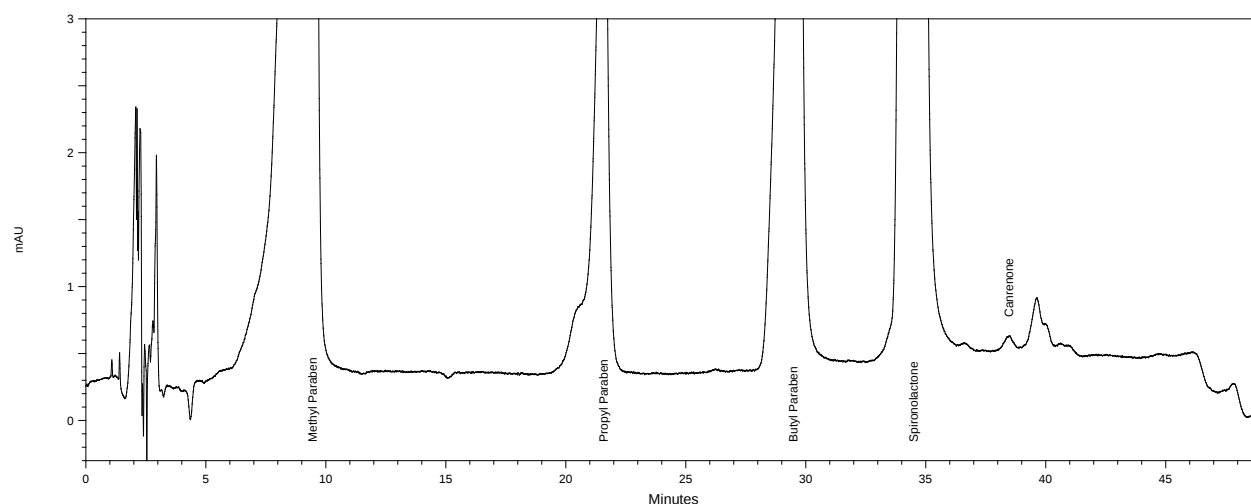

The chromatogram shows four large peaks methyl paraben (off scale at 3 mAU) at 9.0 minutes, propyl paraben (off scale at 3 mAU) at 21.5 minutes and butyl paraben (off scale at 3 mAU) at 29.0 minutes, spironolactone (off scale at 3 mAU) at 34.5 minutes and one small peak canrenone (>1 mAU) at 38.4 minutes.